




## WEIGHING TABLE WITH LEFT GRANITE PLATE

**Length** 800, 900, 1000

**Width** 600, 700, 800

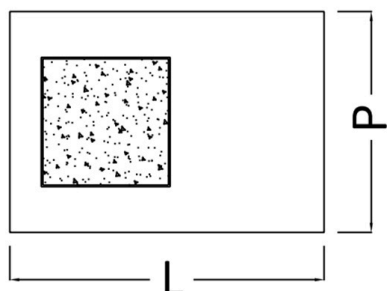


Frame structure in stainless steel AISI 304 fully welded 40x40mm tube, 1.00mm thick, rear legs 80mm advanced, suitable for rounded floor corner. Rear and side welded reinforcement rails. Working top in AISI 304 sheet, 1.50mm thick, with inner plywood reinforcement and stainless steel lower cover. Opening/s on the work surface to allow the housing of 1 or 2 400x400x40Hmm black granite plates, supported by anti-vibration rubber spacers. Rounded (R=10mm) front upper edge, as well as rounded (R=2mm) side and back edges. On adjustable feet.

---

**SKU:** N/A



		W		
		600	700	800
L	Cod.	Cod.	Cod.	Cod.

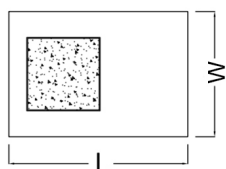

800	T07201001	T07202001	T07203001
900	T07201002	T07202002	T07203002
1000	T07201003	T07202003	T07203003

**W** = Width / Profondità / Breite / Profondeur.

**L** = Length / Larghezza / Länge / Longueur.