




## WEIGHING TABLE WITH REAR UPSTAND WITH DOUBLE GRANITE PLATES

**Length** 1400, 1600, 1800

**Width** 600, 700, 800

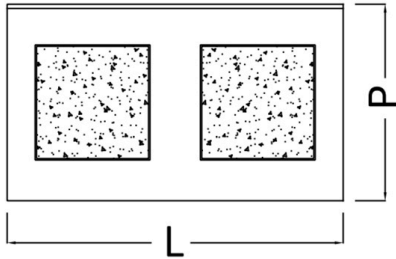

Frame structure in stainless steel AISI 304 fully welded 40x40mm tube, 1.00mm thick, rear legs 80mm advanced, suitable for rounded floor corner. Rear and side welded reinforcement rails. Working top in AISI 304 sheet, 1.50mm thick, with inner plywood reinforcement and stainless steel lower cover. Opening/s on the work surface to allow the housing of 1 or 2 400x400x40Hmm black granite plates, supported by anti-vibration rubber spacers. Rounded (R=10mm) front upper edge, as well as rounded (R=2mm) side and back edges. With rear stainless steel upstand H=100mm, welded and closed on the back, and with R=10mm rounded lower edge. On adjustable feet.

---

**SKU:** N/A



	W		
	600	700	800
L	Cod.	Cod.	Cod.

	1400	T07801001	T07802001	T07803001
	1600	T07801002	T07802002	T07803002
	1800	T07801003	T07802003	T07803003

**W** = Width / Profondità / Breite / Profondeur.

**L** = Length / Larghezza / Länge / Longueur.