




WEIGHING TABLE WITH REAR UPSTAND WITH RIGHT GRANITE PLATE

Length 800, 900, 1000

Width 600, 700, 800

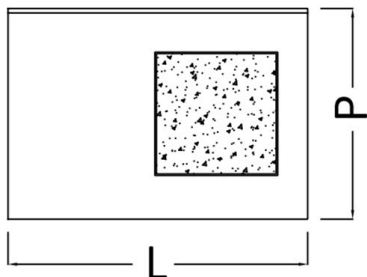


Frame structure in stainless steel AISI 304 fully welded 40x40mm tube, 1.00mm thick, rear legs 80mm advanced, suitable for rounded floor corner. Rear and side welded reinforcement rails. Working top in AISI 304 sheet, 1.50mm thick, with inner plywood reinforcement and stainless steel lower cover. Opening/s on the work surface to allow the housing of 1 or 2 400x400x40Hmm black granite plates, supported by anti-vibration rubber spacers. Rounded (R=10mm) front upper edge, as well as rounded (R=2mm) side and back edges. With rear stainless steel upstand H=100mm, welded and closed on the back, and with R=10mm rounded lower edge. On adjustable feet.

SKU: N/A



| | | W | | |
|---|--|------|------|------|
| | | 600 | 700 | 800 |
| L | | Cod. | Cod. | Cod. |

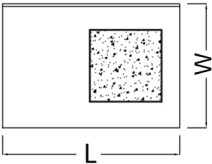

| | | | |
|------|-----------|-----------|-----------|
| 800 | T07601001 | T07602001 | T07603001 |
| 900 | T07601002 | T07602002 | T07603002 |
| 1000 | T07601003 | T07602003 | T07603003 |

W = Width / Profondità / Breite / Profondeur.
L = Length / Larghezza / Länge / Longueur.